



Axxis Segmented

V-TRAK

The versatile solution for complex seating issues





Axxis Segmented

V-TRAK

PERFORMANCE
HEALTH PRODUCTS

The versatile solution for complex seating issues

AXXIS SEGMENTED BACKRESTS - The complexity and variety of the human form frequently requires a wide range of adjustments for both comfort and functional seating.

PAEDIATRIC AND ADULT SIZES FROM ONE RANGE - Axxis segmented backrests allow an exceptional range of paediatric and adult sizes and shapes to be constructed from a range of modular parts. Each segment has all the special features of the Axxis backrests. They have integral lateral supports for increased stability, internal straps for lumbar profiling and for adjusting the resilience and feel of the backrest.

PERMUTATIONS - The modules can be arranged in any combination of size and orientation to suit different body shapes, spine profiles and pressure distribution. Each section can be mounted individually or joined to others. When mounted individually, each segment can be adjusted independently making the system suitable for many complex conditions including Scoliosis.

PRODUCT SELECTION - Quick and convenient selection using the unique V-Trak sizing tape to ensure the correct choice of height, width and combination of backrest segments.

FULL CONTROL - Axxis segmented backrests incorporate all the advantages of the Axxis backrests with complete control of comfort, shape, size, position and resilience.

V-TRAK HARDWARE - At the heart of a good seating system is the ability to provide safe and versatile positioning. With its extensive adjustment in all planes V-Trak provides the ultimate in positioning flexibility.



COMPATIBILITY - Axxis segmented backrests use the same unique mounting and adjustment hardware as Axxis one-piece backrests. This ensures that segmented backrests can be fitted to virtually any wheelchair in common use and will accept all V-Trak accessories and most third party accessories. This gives the client and the clinician the widest possible choice and the opportunity to provide the very best solution.

SAFETY - Axxis backrests mounted on V-Trak hardware are the first backrest systems to pass the new ANSI RESNA Crash Testing Standard WC20 and to be accredited WC20 Compliant.


Axxis Segmented Combinations



Spine profiling
Anterior Pelvic Tilt

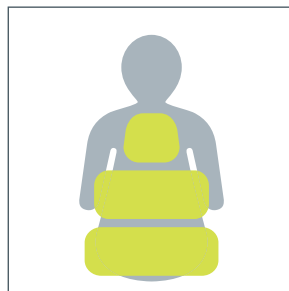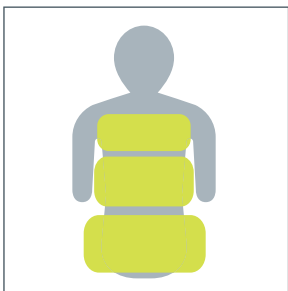


Accommodates complex profiles - **Trunk Taper**

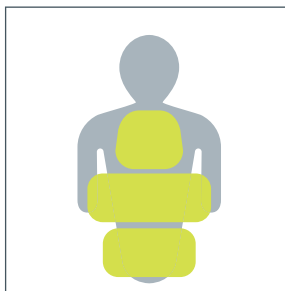


Accommodates complex profiles - **Scoliosis**


Axxis Complex Solutions

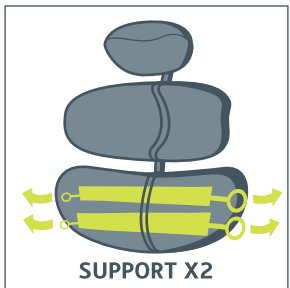


The 2 internal straps of each segment provide adjustable lumbar profile

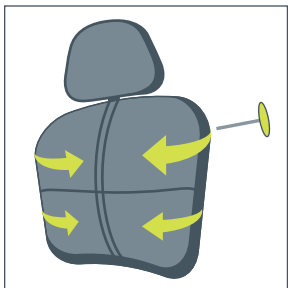

Axxis segments can be combined with our single piece backrests

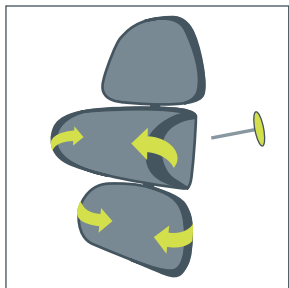

Independently adjustable integral wings for symmetrical and asymmetrical fit

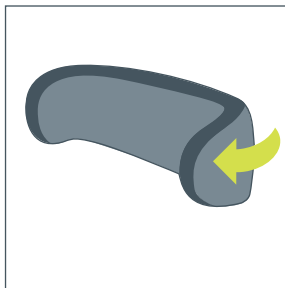

Deep Contouring

Segmented Conditions/Features



Severe Scoliosis



Lumbar Lordosis
Anterior Pelvic Tilt



Pelvic Rotation



Spinal Hernia



Adjustment For
Growth



All Types of
Wheelchair



Crash tested
WC20-2007



Literature &
Training DVD



Easy product
selection

Authorised Distributor:



PERFORMANCE HEALTH PRODUCTS LTD

UNIT 1, WESTSIDE CAMBRIAN INDUSTRIAL ESTATE, COEDCAE LANE, PONTYCLUN, WALES CF72 9EX
TEL +44 (0)1443 236 530 | FAX +44 (0)1443 239 355 | INFO@V-TRAK.COM | WWW.V-TRAK.COM